

```
In [1]: import requests
import pandas as pd
import matplotlib.pyplot as plt

# this is the server which hosts all the patient data.
# in the real-world this would be an EHR vendor
# this is a open-source server with synthetic data for developers to learn a
FHIR_BASE = "http://hapi.fhir.org/baseR4"
```

```
In [2]: # Healthcare applications rarely store all their own data. Instead, they
# exchange information with electronic health record (EHR) systems using
# standards. FHIR (Fast Healthcare Interoperability Resources) is a standard
# based on RESTful web API ideas and allows for retrieving patient informati
# from an EHR AND adding information to an EHR (provided the application has
# authorization)

# In this tutorial, we will walk through how to interact with a public FHIR
# server to get information about a patient.

# We will also create a synthetic patient with clinical information (because
# public servers often don't have the test patients we might need)

# Pretend you are building a new health application for patients which will
# them some health trends and their current medications. The app should be a
# to:
#   - Retrieve health data such as blood pressure, heart rate, and lab resul
#   - Retrieve active medications

# Pull this information back out and display it as tables and charts so pati
# can see trends in their health over time

# By the end of this exercise, you'll have a toy prototype that demonstrates
# a modern health application can connect to a standards-based health data
# source.
```

```
In [3]: # let's query the FHIR server to get information about a patient.
# we'll find up to 5 patients with last name "Smith" (chosen because these
# public servers often use this name)

# Recall that everything in FHIR is a "Resource" (patients, observations, et

resp = requests.get(f"{FHIR_BASE}/Patient", params={"name": "Smith", "_count
patients = resp.json()

for entry in patients.get("entry", []):
    patient = entry["resource"]
    print(patient["id"], patient.get("name"))
```

```
597226 [{'use': 'official', 'family': 'Smith', 'given': ['Sean']}]
597372 [{'family': 'Joe', 'given': ['Smith']}]
597502 [{'family': 'Joe', 'given': ['Smith']}]
597509 [{'family': 'Joe', 'given': ['Smith']}]
618626 [{'use': 'official', 'text': 'Julia Smith', 'family': 'Smith', 'given': ['Julia']}]
```

```
In [4]: # we can see that there are 3 patients with this last name. each has a unique ID
# let's look at "Observations" from the first, limiting to the first 50

patient_id = "pat2" # replace with real ID from above
resp = requests.get(f"{FHIR_BASE}/Observation", params={"patient": patient_id})
obs = resp.json()

# and now let's print out all the observations
obs

# it's a lot! copy into a JSON viewer to look at the structure
# https://jsonviewer.stack.hu/
```

```

Out[4]: {'resourceType': 'Bundle',
        'id': '5ce8ae1d-b358-4f20-813c-710fe72ec8f0',
        'meta': {'lastUpdated': '2025-10-14T19:32:25.539+00:00'},
        'type': 'searchset',
        'link': [{'relation': 'self',
                  'url': 'https://hapi.fhir.org/baseR4/0bservation?_count=50&patient=pat2'}],
        {'relation': 'next',
         'url': 'https://hapi.fhir.org/baseR4?_getpages=5ce8ae1d-b358-4f20-813c-710fe72ec8f0&_getpagesoffset=50&_count=50&_pretty=true&_bundletype=searchset'}],
        'entry': [{'fullUrl': 'https://hapi.fhir.org/baseR4/0bservation/697307',
                    'resource': {'resourceType': 'Observation',
                                'id': '697307',
                                'meta': {'versionId': '1',
                                          'lastUpdated': '2020-03-21T13:35:12.175+00:00',
                                          'source': '#jRCIxW30Q18JIL0r'},
                                'identifier': [{'system': 'http://acme.com/lab/observation',
                                              'value': 'ObserationId_5234342'}],
                                'status': 'final',
                                'code': {'coding': [{'system': 'http://loinc.org',
                                                    'code': '35200-5',
                                                    'display': 'Cholesterol [Moles/\u200bvolume] in Serum or Plasma'}],
                                          'text': 'Cholesterol'},
                                'subject': {'reference': 'Patient/pat2'},
                                'performer': [{'reference': 'Practitioner/697300'}],
                                'valueQuantity': {'value': 6.3,
                                                  'unit': 'mmol/L',
                                                  'system': 'http://unitsofmeasure.org',
                                                  'code': 'mmol/L'},
                                'referenceRange': [{'high': {'value': 1.5,
                                                            'unit': 'mmol/L',
                                                            'system': 'http://unitsofmeasure.org',
                                                            'code': 'mmol/L'}}]},
                    'search': {'mode': 'match'}}],
        {'fullUrl': 'https://hapi.fhir.org/baseR4/0bservation/2138112',
         'resource': {'resourceType': 'Observation',
                     'id': '2138112',
                     'meta': {'versionId': '1',
                               'lastUpdated': '2021-05-25T12:52:24.320+00:00',
                               'source': '#RSHqrS7oI1Tx51iU'},
                     'status': 'final',
                     'code': {'coding': [{'system': 'http://loinc.org',
                                           'code': '718-7',
                                           'display': 'Hemoglobin [Mass/volume] in Blood'}],
                               'text': 'Haemoglobin'},
                     'subject': {'reference': 'Patient/pat2'},
                     'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9879522f',
                                    'display': 'Acme Laboratory, Inc'}],
                     'valueQuantity': {'value': 176,
                                       'unit': 'g/L',
                                       'system': 'http://unitsofmeasure.org',
                                       'code': 'g/L'},
                     'referenceRange': [{'low': {'value': 135,
                                                'unit': 'g/L',

```

```

      'system': 'http://unitsofmeasure.org',
      'code': 'g/L'},
    'high': {'value': 180,
      'unit': 'g/L',
      'system': 'http://unitsofmeasure.org',
      'code': 'g/L'}}}],
    'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138113',
      'resource': {'resourceType': 'Observation',
        'id': '2138113',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:52:24.320+00:00',
          'source': '#RSHqrS7oI1Tx51iU'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '789-8',
          'display': 'Erythrocytes [# /volume] in Blood by Automated count'}]},
        'text': 'Red Cell Count'},
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}]},
        'valueQuantity': {'value': 5.9,
          'unit': 'x10*12/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*12/L'},
        'referenceRange': [{'low': {'value': 4.2,
          'unit': 'x10*12/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*12/L'},
          'high': {'value': 6.0,
          'unit': 'x10*12/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*12/L'}}]}],
        'search': {'mode': 'match'}},
        {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138114',
          'resource': {'resourceType': 'Observation',
            'id': '2138114',
            'meta': {'versionId': '1',
              'lastUpdated': '2021-05-25T12:52:24.320+00:00',
              'source': '#RSHqrS7oI1Tx51iU'},
            'status': 'final',
            'code': {'coding': [{'system': 'http://loinc.org',
              'code': '4544-3',
              'display': 'Hematocrit [Volume Fraction] of Blood by Automated coun
t'}]},
            'text': 'Haematocrit'},
            'subject': {'reference': 'Patient/pat2'},
            'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
              'display': 'Acme Laboratory, Inc'}]},
            'valueQuantity': {'value': 55, 'unit': '%'},
            'interpretation': [{'coding': [{'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
              'code': 'H'}}]}],
            'referenceRange': [{'low': {'value': 38, 'unit': '%'},

```

```
      'high': {'value': 52, 'unit': '%'}}}},
    'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138115',
      'resource': {'resourceType': 'Observation',
        'id': '2138115',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:52:24.320+00:00',
          'source': '#RSHqrS7oI1Tx51iU'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '787-2',
          'display': 'Erythrocyte mean corpuscular volume [Entitic volume] by
Automated count'}]},
        'text': 'Mean Cell Volume',
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}]},
        'valueQuantity': {'value': 99,
          'unit': 'fL',
          'system': 'http://unitsofmeasure.org',
          'code': 'fL'},
        'interpretation': [{'coding': [{'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
          'code': 'H'}}]}]},
        'referenceRange': [{'low': {'value': 80,
          'unit': 'fL',
          'system': 'http://unitsofmeasure.org',
          'code': 'fL'},
          'high': {'value': 98,
            'unit': 'fL',
            'system': 'http://unitsofmeasure.org',
            'code': 'fL'}}]}]},
    'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138116',
      'resource': {'resourceType': 'Observation',
        'id': '2138116',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:52:24.320+00:00',
          'source': '#RSHqrS7oI1Tx51iU'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '785-6',
          'display': 'Erythrocyte mean corpuscular hemoglobin [Entitic mass] b
y Automated count'}]},
        'text': 'Mean Cell Haemoglobin',
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}]},
        'valueQuantity': {'value': 36,
          'unit': 'pg',
          'system': 'http://unitsofmeasure.org',
          'code': 'pg'},
        'interpretation': [{'coding': [{'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
```

```
      'code': 'H'}}]]],
    'referenceRange': [{ 'low': { 'value': 27,
      'unit': 'pg',
      'system': 'http://unitsofmeasure.org',
      'code': 'pg'},
    'high': { 'value': 35,
      'unit': 'pg',
      'system': 'http://unitsofmeasure.org',
      'code': 'pg'}}]],
    'search': { 'mode': 'match' }},
{ 'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138117',
  'resource': { 'resourceType': 'Observation',
    'id': '2138117',
    'meta': { 'versionId': '1',
      'lastUpdated': '2021-05-25T12:52:24.320+00:00',
      'source': '#RSHqrS7oI1Tx51iU'},
    'status': 'final',
    'code': { 'coding': [{ 'system': 'http://loinc.org',
      'code': '777-3',
      'display': 'Platelets [# /volume] in Blood by Automated count'}],
      'text': 'Platelet Count'},
    'subject': { 'reference': 'Patient/pat2'},
    'performer': [{ 'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
      'display': 'Acme Laboratory, Inc'}],
    'valueQuantity': { 'value': 444,
      'unit': 'x10*9/L',
      'system': 'http://unitsofmeasure.org',
      'code': '10*9/L'},
    'referenceRange': [{ 'low': { 'value': 150,
      'unit': 'x10*9/L',
      'system': 'http://unitsofmeasure.org',
      'code': '10*9/L'},
    'high': { 'value': 450,
      'unit': 'x10*9/L',
      'system': 'http://unitsofmeasure.org',
      'code': '10*9/L'}}]],
    'search': { 'mode': 'match' }},
{ 'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138118',
  'resource': { 'resourceType': 'Observation',
    'id': '2138118',
    'meta': { 'versionId': '1',
      'lastUpdated': '2021-05-25T12:52:24.320+00:00',
      'source': '#RSHqrS7oI1Tx51iU'},
    'status': 'final',
    'code': { 'coding': [{ 'system': 'http://loinc.org',
      'code': '6690-2',
      'display': 'Leukocytes [# /volume] in Blood by Automated count'}],
      'text': 'White Cell Count'},
    'subject': { 'reference': 'Patient/pat2'},
    'performer': [{ 'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
      'display': 'Acme Laboratory, Inc'}],
    'valueQuantity': { 'value': 4.6,
      'unit': 'x10*9/L',
      'system': 'http://unitsofmeasure.org',
```

```

      'code': '10*9/L'},
      'referenceRange': [{ 'low': { 'value': 4.0,
        'unit': 'x10*9/L',
        'system': 'http://unitsofmeasure.org',
        'code': '10*9/L'},
        'high': { 'value': 11.0,
        'unit': 'x10*9/L',
        'system': 'http://unitsofmeasure.org',
        'code': '10*9/L'}}}],
      'search': { 'mode': 'match' }},
    { 'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138119',
      'resource': { 'resourceType': 'Observation',
        'id': '2138119',
        'meta': { 'versionId': '1',
          'lastUpdated': '2021-05-25T12:52:24.320+00:00',
          'source': '#RSHqrS7oI1Tx51iU'},
        'status': 'final',
        'code': { 'coding': [{ 'system': 'http://loinc.org',
          'code': '770-8',
          'display': 'Neutrophils/100 leukocytes in Blood by Automated coun
t'}],
          'text': 'Neutrophils'},
        'subject': { 'reference': 'Patient/pat2'},
        'performer': [{ 'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}],
        'valueQuantity': { 'value': 20,
          'unit': '%',
          'system': 'http://unitsofmeasure.org',
          'code': '%'}},
        'search': { 'mode': 'match' }},
    { 'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138120',
      'resource': { 'resourceType': 'Observation',
        'id': '2138120',
        'meta': { 'versionId': '1',
          'lastUpdated': '2021-05-25T12:52:24.320+00:00',
          'source': '#RSHqrS7oI1Tx51iU'},
        'status': 'final',
        'code': { 'coding': [{ 'system': 'http://loinc.org',
          'code': '751-8',
          'display': 'Neutrophils [# /volume] in Blood by Automated count'}],
          'text': 'Neutrophils'},
        'subject': { 'reference': 'Patient/pat2'},
        'performer': [{ 'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}],
        'valueQuantity': { 'value': 0.9,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L'},
        'interpretation': [{ 'coding': [{ 'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
          'code': 'LL'}]}],
        'referenceRange': [{ 'low': { 'value': 2.0,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',

```

```
      'code': '10*9/L'},
      'high': {'value': 7.5,
        'unit': 'x10*9/L',
        'system': 'http://unitsofmeasure.org',
        'code': '10*9/L'}}}],
    'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138121',
      'resource': {'resourceType': 'Observation',
        'id': '2138121',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:52:24.320+00:00',
          'source': '#RSHqrS7oI1Tx51iU'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '736-9',
          'display': 'Lymphocytes/100 leukocytes in Blood by Automated coun
t'}]},
        'text': 'Lymphocytes'},
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}]},
        'valueQuantity': {'value': 20,
          'unit': '%',
          'system': 'http://unitsofmeasure.org',
          'code': '%'}},
        'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138122',
      'resource': {'resourceType': 'Observation',
        'id': '2138122',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:52:24.320+00:00',
          'source': '#RSHqrS7oI1Tx51iU'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '731-0',
          'display': 'Lymphocytes [# /volume] in Blood by Automated count'}]},
        'text': 'Lymphocytes'},
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}]},
        'valueQuantity': {'value': 0.9,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L'},
        'interpretation': [{'coding': [{'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
          'code': 'L'}]}]},
        'referenceRange': [{'low': {'value': 1.1,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L'},
        'high': {'value': 4.0,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
```



```

      'code': '10*9/L'}}}],
    'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138123',
      'resource': {'resourceType': 'Observation',
        'id': '2138123',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:52:24.320+00:00',
          'source': '#RSHqrS7oI1Tx51iU'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '5905-5',
          'display': 'Monocytes/100 leukocytes in Blood by Automated count'}]},
        'text': 'Monocytes'},
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}]},
        'valueQuantity': {'value': 20,
          'unit': '%',
          'system': 'http://unitsofmeasure.org',
          'code': '%'}},
        'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138124',
      'resource': {'resourceType': 'Observation',
        'id': '2138124',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:52:24.320+00:00',
          'source': '#RSHqrS7oI1Tx51iU'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '742-7',
          'display': 'Monocytes [# /volume] in Blood by Automated count'}]},
        'text': 'Monocytes'},
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}]},
        'valueQuantity': {'value': 0.9,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L'},
        'referenceRange': [{'low': {'value': 0.2,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L'},
            'high': {'value': 1.0,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L'}}}],
        'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138125',
      'resource': {'resourceType': 'Observation',
        'id': '2138125',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:52:24.320+00:00',
          'source': '#RSHqrS7oI1Tx51iU'},

```

```
'status': 'final',
'code': {'coding': [{'system': 'http://loinc.org',
  'code': '713-8',
  'display': 'Eosinophils/100 leukocytes in Blood by Automated coun
t'}]},
'text': 'Eosinophils'},
'subject': {'reference': 'Patient/pat2'},
'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
  'display': 'Acme Laboratory, Inc'}]},
'valueQuantity': {'value': 20,
  'unit': '%',
  'system': 'http://unitsofmeasure.org',
  'code': '%'}},
'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138126',
'resource': {'resourceType': 'Observation',
'id': '2138126',
'meta': {'versionId': '1',
  'lastUpdated': '2021-05-25T12:52:24.320+00:00',
  'source': '#RSHqrS7oI1Tx51iU'},
'status': 'final',
'code': {'coding': [{'system': 'http://loinc.org',
  'code': '711-2',
  'display': 'Eosinophils [# /volume] in Blood by Automated count'}]},
'text': 'Eosinophils'},
'subject': {'reference': 'Patient/pat2'},
'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
  'display': 'Acme Laboratory, Inc'}]},
'valueQuantity': {'value': 0.92,
  'unit': 'x10*9/L',
  'system': 'http://unitsofmeasure.org',
  'code': '10*9/L'},
'interpretation': [{'coding': [{'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
  'code': 'HH'}]}]},
'referenceRange': [{'low': {'value': 0.04,
  'unit': 'x10*9/L',
  'system': 'http://unitsofmeasure.org',
  'code': '10*9/L'},
  'high': {'value': 0.4,
  'unit': 'x10*9/L',
  'system': 'http://unitsofmeasure.org',
  'code': '10*9/L'}}]},
'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138127',
'resource': {'resourceType': 'Observation',
'id': '2138127',
'meta': {'versionId': '1',
  'lastUpdated': '2021-05-25T12:52:24.320+00:00',
  'source': '#RSHqrS7oI1Tx51iU'},
'status': 'final',
'code': {'coding': [{'system': 'http://loinc.org',
  'code': '706-2',
  'display': 'Basophils/100 leukocytes in Blood by Automated count'}]},
```

```
      'text': 'Basophils'},
      'subject': {'reference': 'Patient/pat2'},
      'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
        'display': 'Acme Laboratory, Inc'}],
      'valueQuantity': {'value': 20,
        'unit': '%',
        'system': 'http://unitsofmeasure.org',
        'code': '%'}},
      'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138128',
      'resource': {'resourceType': 'Observation',
        'id': '2138128',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:52:24.320+00:00',
          'source': '#RSHqrS7oI1Tx51iU'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '704-7',
          'display': 'Basophils [# /volume] in Blood by Automated count'}]},
        'text': 'Basophils'},
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}],
        'valueQuantity': {'value': 0.92,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L'},
        'referenceRange': [{'high': {'value': 0.21,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L'}}]},
        'search': {'mode': 'match'}},
      'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138324',
      'resource': {'resourceType': 'Observation',
        'id': '2138324',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:58:32.061+00:00',
          'source': '#CZuTWAJ8dajUPU6v'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '718-7',
          'display': 'Emoglobina'}]},
        'text': 'Emoglobina'},
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}],
        'valueQuantity': {'value': 176,
          'unit': 'g/L',
          'system': 'http://unitsofmeasure.org',
          'code': 'g/L'},
        'referenceRange': [{'low': {'value': 135,
          'unit': 'g/L',
          'system': 'http://unitsofmeasure.org',
```

```

      'code': 'g/L'},
      'high': {'value': 180,
        'unit': 'g/L',
        'system': 'http://unitsofmeasure.org',
        'code': 'g/L'}}}],
    'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138325',
      'resource': {'resourceType': 'Observation',
        'id': '2138325',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:58:32.061+00:00',
          'source': '#CZuTWAJ8dajUPU6v'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '789-8',
          'display': 'Erythrocytes [#/volume] in Blood by Automated count'}]},
        'text': 'Red Cell Count',
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}]},
        'valueQuantity': {'value': 5.9,
          'unit': 'x10*12/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*12/L'},
        'referenceRange': [{'low': {'value': 4.2,
          'unit': 'x10*12/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*12/L'},
          'high': {'value': 6.0,
            'unit': 'x10*12/L',
            'system': 'http://unitsofmeasure.org',
            'code': '10*12/L'}}]}],
        'search': {'mode': 'match'}},
      {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138326',
        'resource': {'resourceType': 'Observation',
          'id': '2138326',
          'meta': {'versionId': '1',
            'lastUpdated': '2021-05-25T12:58:32.061+00:00',
            'source': '#CZuTWAJ8dajUPU6v'},
          'status': 'final',
          'code': {'coding': [{'system': 'http://loinc.org',
            'code': '4544-3',
            'display': 'Hematocrit [Volume Fraction] of Blood by Automated coun
t'}]},
          'text': 'Haematocrit',
          'subject': {'reference': 'Patient/pat2'},
          'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
            'display': 'Acme Laboratory, Inc'}]},
          'valueQuantity': {'value': 55, 'unit': '%'},
          'interpretation': [{'coding': [{'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
            'code': 'H'}}]}],
          'referenceRange': [{'low': {'value': 38, 'unit': '%'},
            'high': {'value': 52, 'unit': '%'}}]}],

```

```
'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138327',
 'resource': {'resourceType': 'Observation',
  'id': '2138327',
  'meta': {'versionId': '1',
   'lastUpdated': '2021-05-25T12:58:32.061+00:00',
   'source': '#CZuTWAJ8dajUPU6v'},
  'status': 'final',
  'code': {'coding': [{'system': 'http://loinc.org',
   'code': '787-2',
   'display': 'Erythrocyte mean corpuscular volume [Entitic volume] by
Automated count'}]},
  'text': 'Mean Cell Volume',
  'subject': {'reference': 'Patient/pat2'},
  'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
   'display': 'Acme Laboratory, Inc'}]},
  'valueQuantity': {'value': 99,
   'unit': 'fL',
   'system': 'http://unitsofmeasure.org',
   'code': 'fL'},
  'interpretation': [{'coding': [{'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
   'code': 'H'}]}]},
  'referenceRange': [{'low': {'value': 80,
   'unit': 'fL',
   'system': 'http://unitsofmeasure.org',
   'code': 'fL'},
   'high': {'value': 98,
   'unit': 'fL',
   'system': 'http://unitsofmeasure.org',
   'code': 'fL'}}]}],
'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138328',
 'resource': {'resourceType': 'Observation',
  'id': '2138328',
  'meta': {'versionId': '1',
   'lastUpdated': '2021-05-25T12:58:32.061+00:00',
   'source': '#CZuTWAJ8dajUPU6v'},
  'status': 'final',
  'code': {'coding': [{'system': 'http://loinc.org',
   'code': '785-6',
   'display': 'Erythrocyte mean corpuscular hemoglobin [Entitic mass] b
y Automated count'}]},
  'text': 'Mean Cell Haemoglobin',
  'subject': {'reference': 'Patient/pat2'},
  'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
   'display': 'Acme Laboratory, Inc'}]},
  'valueQuantity': {'value': 36,
   'unit': 'pg',
   'system': 'http://unitsofmeasure.org',
   'code': 'pg'},
  'interpretation': [{'coding': [{'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
   'code': 'H'}]}]},
```

```
'referenceRange': [{ 'low': { 'value': 27,
  'unit': 'pg',
  'system': 'http://unitsofmeasure.org',
  'code': 'pg' },
  'high': { 'value': 35,
  'unit': 'pg',
  'system': 'http://unitsofmeasure.org',
  'code': 'pg' } }],
'search': { 'mode': 'match' },
{ 'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138329',
  'resource': { 'resourceType': 'Observation',
    'id': '2138329',
    'meta': { 'versionId': '1',
      'lastUpdated': '2021-05-25T12:58:32.061+00:00',
      'source': '#CZuTWAJ8dajUPU6v' },
    'status': 'final',
    'code': { 'coding': [{ 'system': 'http://loinc.org',
      'code': '777-3',
      'display': 'Platelets [# /volume] in Blood by Automated count' } ] },
    'text': 'Platelet Count',
    'subject': { 'reference': 'Patient/pat2' },
    'performer': [{ 'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
      'display': 'Acme Laboratory, Inc' } ] },
    'valueQuantity': { 'value': 444,
      'unit': 'x10*9/L',
      'system': 'http://unitsofmeasure.org',
      'code': '10*9/L' },
    'referenceRange': [{ 'low': { 'value': 150,
      'unit': 'x10*9/L',
      'system': 'http://unitsofmeasure.org',
      'code': '10*9/L' },
      'high': { 'value': 450,
      'unit': 'x10*9/L',
      'system': 'http://unitsofmeasure.org',
      'code': '10*9/L' } } } ],
    'search': { 'mode': 'match' } },
{ 'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138330',
  'resource': { 'resourceType': 'Observation',
    'id': '2138330',
    'meta': { 'versionId': '1',
      'lastUpdated': '2021-05-25T12:58:32.061+00:00',
      'source': '#CZuTWAJ8dajUPU6v' },
    'status': 'final',
    'code': { 'coding': [{ 'system': 'http://loinc.org',
      'code': '6690-2',
      'display': 'Leukocytes [# /volume] in Blood by Automated count' } ] },
    'text': 'White Cell Count',
    'subject': { 'reference': 'Patient/pat2' },
    'performer': [{ 'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
      'display': 'Acme Laboratory, Inc' } ] },
    'valueQuantity': { 'value': 4.6,
      'unit': 'x10*9/L',
      'system': 'http://unitsofmeasure.org',
      'code': '10*9/L' },
```

```

      'referenceRange': [{ 'low': { 'value': 4.0,
        'unit': 'x10*9/L',
        'system': 'http://unitsofmeasure.org',
        'code': '10*9/L' },
        'high': { 'value': 11.0,
        'unit': 'x10*9/L',
        'system': 'http://unitsofmeasure.org',
        'code': '10*9/L' } } ],
      'search': { 'mode': 'match' } },
    { 'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138331',
      'resource': { 'resourceType': 'Observation',
        'id': '2138331',
        'meta': { 'versionId': '1',
          'lastUpdated': '2021-05-25T12:58:32.061+00:00',
          'source': '#CZuTWAJ8dajUPU6v' },
        'status': 'final',
        'code': { 'coding': [{ 'system': 'http://loinc.org',
          'code': '770-8',
          'display': 'Neutrophils/100 leukocytes in Blood by Automated coun
t' } ] },
        'text': 'Neutrophils' },
        'subject': { 'reference': 'Patient/pat2' },
        'performer': [{ 'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc' } ] },
        'valueQuantity': { 'value': 20,
          'unit': '%',
          'system': 'http://unitsofmeasure.org',
          'code': '%'} },
        'search': { 'mode': 'match' } },
    { 'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138332',
      'resource': { 'resourceType': 'Observation',
        'id': '2138332',
        'meta': { 'versionId': '1',
          'lastUpdated': '2021-05-25T12:58:32.061+00:00',
          'source': '#CZuTWAJ8dajUPU6v' },
        'status': 'final',
        'code': { 'coding': [{ 'system': 'http://loinc.org',
          'code': '751-8',
          'display': 'Neutrophils [# /volume] in Blood by Automated count' } ] },
        'text': 'Neutrophils' },
        'subject': { 'reference': 'Patient/pat2' },
        'performer': [{ 'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc' } ] },
        'valueQuantity': { 'value': 0.9,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L' },
        'interpretation': [{ 'coding': [{ 'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
          'code': 'LL' } } ] },
        'referenceRange': [{ 'low': { 'value': 2.0,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L' },

```

```

      'high': {'value': 7.5,
        'unit': 'x10*9/L',
        'system': 'http://unitsofmeasure.org',
        'code': '10*9/L'}}}],
    'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138333',
      'resource': {'resourceType': 'Observation',
        'id': '2138333',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:58:32.061+00:00',
          'source': '#CZuTWAJ8dajUPU6v'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '736-9',
          'display': 'Lymphocytes/100 leukocytes in Blood by Automated coun
t'}]},
        'text': 'Lymphocytes'},
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}]},
        'valueQuantity': {'value': 20,
          'unit': '%',
          'system': 'http://unitsofmeasure.org',
          'code': '%'}},
        'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138334',
      'resource': {'resourceType': 'Observation',
        'id': '2138334',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:58:32.061+00:00',
          'source': '#CZuTWAJ8dajUPU6v'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '731-0',
          'display': 'Lymphocytes [# /volume] in Blood by Automated count'}]},
        'text': 'Lymphocytes'},
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}]},
        'valueQuantity': {'value': 0.9,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L'},
        'interpretation': [{'coding': [{'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
          'code': 'L'}]}]},
        'referenceRange': [{'low': {'value': 1.1,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L'},
          'high': {'value': 4.0,
            'unit': 'x10*9/L',
            'system': 'http://unitsofmeasure.org',
            'code': '10*9/L'}}}],

```



```
'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138335',
 'resource': {'resourceType': 'Observation',
  'id': '2138335',
  'meta': {'versionId': '1',
   'lastUpdated': '2021-05-25T12:58:32.061+00:00',
   'source': '#CZuTWAJ8dajUPU6v'},
  'status': 'final',
  'code': {'coding': [{'system': 'http://loinc.org',
   'code': '5905-5',
   'display': 'Monocytes/100 leukocytes in Blood by Automated count'}]},
  'text': 'Monocytes'},
  'subject': {'reference': 'Patient/pat2'},
  'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
   'display': 'Acme Laboratory, Inc'}]},
  'valueQuantity': {'value': 20,
   'unit': '%',
   'system': 'http://unitsofmeasure.org',
   'code': '%'}},
 'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138336',
 'resource': {'resourceType': 'Observation',
  'id': '2138336',
  'meta': {'versionId': '1',
   'lastUpdated': '2021-05-25T12:58:32.061+00:00',
   'source': '#CZuTWAJ8dajUPU6v'},
  'status': 'final',
  'code': {'coding': [{'system': 'http://loinc.org',
   'code': '742-7',
   'display': 'Monocytes [# /volume] in Blood by Automated count'}]},
  'text': 'Monocytes'},
  'subject': {'reference': 'Patient/pat2'},
  'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
   'display': 'Acme Laboratory, Inc'}]},
  'valueQuantity': {'value': 0.9,
   'unit': 'x10*9/L',
   'system': 'http://unitsofmeasure.org',
   'code': '10*9/L'},
  'referenceRange': [{'low': {'value': 0.2,
   'unit': 'x10*9/L',
   'system': 'http://unitsofmeasure.org',
   'code': '10*9/L'},
   'high': {'value': 1.0,
   'unit': 'x10*9/L',
   'system': 'http://unitsofmeasure.org',
   'code': '10*9/L'}}]},
 'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138337',
 'resource': {'resourceType': 'Observation',
  'id': '2138337',
  'meta': {'versionId': '1',
   'lastUpdated': '2021-05-25T12:58:32.061+00:00',
   'source': '#CZuTWAJ8dajUPU6v'},
  'status': 'final',
```

```

      'code': {'coding': [{'system': 'http://loinc.org',
        'code': '713-8',
        'display': 'Eosinophils/100 leukocytes in Blood by Automated coun
t'}]},
      'text': 'Eosinophils'},
      'subject': {'reference': 'Patient/pat2'},
      'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
        'display': 'Acme Laboratory, Inc'}]},
      'valueQuantity': {'value': 20,
        'unit': '%',
        'system': 'http://unitsofmeasure.org',
        'code': '%'}},
      'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138338',
      'resource': {'resourceType': 'Observation',
        'id': '2138338',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:58:32.061+00:00',
          'source': '#CZuTWAJ8dajUPU6v'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '711-2',
          'display': 'Eosinophils [# /volume] in Blood by Automated count'}]},
        'text': 'Eosinophils'},
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}]},
        'valueQuantity': {'value': 0.92,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L'},
        'interpretation': [{'coding': [{'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
          'code': 'HH'}]}]},
        'referenceRange': [{'low': {'value': 0.04,
          'unit': 'x10*9/L',
          'system': 'http://unitsofmeasure.org',
          'code': '10*9/L'},
          'high': {'value': 0.4,
            'unit': 'x10*9/L',
            'system': 'http://unitsofmeasure.org',
            'code': '10*9/L'}}]}]},
        'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138339',
      'resource': {'resourceType': 'Observation',
        'id': '2138339',
        'meta': {'versionId': '1',
          'lastUpdated': '2021-05-25T12:58:32.061+00:00',
          'source': '#CZuTWAJ8dajUPU6v'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '706-2',
          'display': 'Basophils/100 leukocytes in Blood by Automated count'}]},
        'text': 'Basophils'},

```

```
'subject': {'reference': 'Patient/pat2'},
'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
  'display': 'Acme Laboratory, Inc'}],
'valueQuantity': {'value': 20,
  'unit': '%',
  'system': 'http://unitsofmeasure.org',
  'code': '%'}},
'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2138340',
  'resource': {'resourceType': 'Observation',
    'id': '2138340',
    'meta': {'versionId': '1',
      'lastUpdated': '2021-05-25T12:58:32.061+00:00',
      'source': '#CZuTWAJ8dajUPU6v'},
    'status': 'final',
    'code': {'coding': [{'system': 'http://loinc.org',
      'code': '704-7',
      'display': 'Basophils [# /volume] in Blood by Automated count'}],
    'text': 'Basophils'},
    'subject': {'reference': 'Patient/pat2'},
    'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
      'display': 'Acme Laboratory, Inc'}],
    'valueQuantity': {'value': 0.92,
      'unit': 'x10*9/L',
      'system': 'http://unitsofmeasure.org',
      'code': '10*9/L'},
    'referenceRange': [{'high': {'value': 0.21,
      'unit': 'x10*9/L',
      'system': 'http://unitsofmeasure.org',
      'code': '10*9/L'}}]},
    'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/2569857',
  'resource': {'resourceType': 'Observation',
    'id': '2569857',
    'meta': {'versionId': '1',
      'lastUpdated': '2021-10-05T07:22:39.267+00:00',
      'source': '#06z0GtJxs0G7L9hv'},
    'subject': {'reference': 'Patient/pat2'}},
    'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801536',
  'resource': {'resourceType': 'Observation',
    'id': '44801536',
    'meta': {'versionId': '1',
      'lastUpdated': '2024-07-11T15:44:08.135+00:00',
      'source': '#tYPNjWj0MUUsuCAKv'},
    'text': {'status': 'generated',
      'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narr
ative: Observation</b><a name="o14"> </a></p><div style="display: inline-bl
ock; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px soli
d #8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0
px">Resource Observation &quot;o14&quot; </p></div><p><b>status</b>: final
</p><p><b>code</b>: Alanine aminotransferase [Enzymatic activity/volume] in
Serum or Plasma <span style="background: LightGoldenRodYellow; margin: 4px;
border: 1px solid khaki"> (<a href="https://loinc.org/">LOINC</a>#1742-6; <
```

```

a href="https://browser.ihtsdotools.org/">SNOMED CT</a>#34608000)</span></p>
<p><b>subject</b>: <a href="patient-example-b.html">Patient/pat2</a> &quot;
;Duck DONALD&quot;</p><p><b>performer</b>: <a href="organization-example-la
b.html">Organization/1832473e-2fe0-452d-abe9-3cdb9879522f: Acme Laboratory,
Inc</a> &quot;Clinical Lab&quot;</p><p><b>value</b>: 54 U/L</p><p><b>interp
retation</b>: High <span style="background: LightGoldenRodYellow; margin: 4
px; border: 1px solid khaki"> (<a href="http://terminology.hl7.org/5.1.0/Co
deSystem-v3-ObservationInterpretation.html">ObservationInterpretation</a>#
H)</span></p><h3>ReferenceRanges</h3><table class="grid"><tr><td>-</td><td>
<b>Low</b></td><td><b>High</b></td></tr><tr><td>*</td><td>0 U/L</td><td>45
U/L</td></tr></table></div>'},
  'status': 'final',
  'code': {'coding': [{'system': 'http://loinc.org',
    'code': '1742-6',
    'display': 'Alanine aminotransferase [Enzymatic activity/volume] in
Serum or Plasma'}],
    {'system': 'http://snomed.info/sct', 'code': '34608000'}}},
  'subject': {'reference': 'Patient/pat2'},
  'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
    'display': 'Acme Laboratory, Inc'}],
  'valueQuantity': {'value': 54, 'unit': 'U/L'},
  'interpretation': [{'coding': [{'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
    'code': 'H'}]}],
  'referenceRange': [{'low': {'value': 0, 'unit': 'U/L'},
    'high': {'value': 45, 'unit': 'U/L'}}],
  'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801537',
  'resource': {'resourceType': 'Observation',
    'id': '44801537',
    'meta': {'versionId': '1',
      'lastUpdated': '2024-07-11T15:44:08.135+00:00',
      'source': '#tYPNjWj0MUuCAKv'},
    'text': {'status': 'generated',
      'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narr
ative: Observation</b><a name="o15"> </a></p><div style="display: inline-bl
ock; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px soli
d #8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0
px">Resource Observation &quot;o15&quot; </p></div><p><b>status</b>: final
</p><p><b>code</b>: Aspartate aminotransferase [Enzymatic activity/volume]
in Serum or Plasma <span style="background: LightGoldenRodYellow; margin: 4
px; border: 1px solid khaki"> (<a href="https://loinc.org/">LOINC</a>#192
0-8; <a href="https://browser.ihtsdotools.org/">SNOMED CT</a>#45896001)</sp
an></p><p><b>subject</b>: <a href="patient-example-b.html">Patient/pat2</a>
&quot;Duck DONALD&quot;</p><p><b>performer</b>: <a href="organization-examp
le-lab.html">Organization/1832473e-2fe0-452d-abe9-3cdb9879522f: Acme Labora
tory, Inc</a> &quot;Clinical Lab&quot;</p><p><b>value</b>: 30 U/L</p><h3>Re
ferenceRanges</h3><table class="grid"><tr><td>-</td><td><b>Low</b></td><td>
<b>High</b></td></tr><tr><td>*</td><td>0 U/L</td><td>41 U/L</td></tr></tabl
e></div>'},
  'status': 'final',
  'code': {'coding': [{'system': 'http://loinc.org',
    'code': '1920-8',
    'display': 'Aspartate aminotransferase [Enzymatic activity/volume] i
n Serum or Plasma'}],

```

```
{'system': 'http://snomed.info/sct', 'code': '45896001'}}]},
'subject': {'reference': 'Patient/pat2'},
'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
'display': 'Acme Laboratory, Inc'}],
'valueQuantity': {'value': 30, 'unit': 'U/L'},
'referenceRange': [{'low': {'value': 0, 'unit': 'U/L'},
'high': {'value': 41, 'unit': 'U/L'}}]},
'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801538',
'resource': {'resourceType': 'Observation',
'id': '44801538',
'meta': {'versionId': '1',
'lastUpdated': '2024-07-11T15:44:08.135+00:00',
'source': '#tYPNjWj0MUUsuCAKv'},
'text': {'status': 'generated',
'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narr
ative: Observation</b><a name="o6"> </a></p><div style="display: inline-blo
ck; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px solid
#8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0px
">Resource Observation &quot;o6&quot; </p></div><p><b>status</b>: final</p>
<p><b>code</b>: Glucose [Moles/volume] in Serum or Plasma <span style="back
ground: LightGoldenRodYellow; margin: 4px; border: 1px solid khaki"> (<a hr
ef="https://loinc.org/">LOINC</a>#14749-6; <a href="https://browser.ihtsdot
ools.org/">SNOMED CT</a>#36048009)</span></p><p><b>subject</b>: <a href="pa
tient-example-b.html">Patient/pat2</a> &quot;Duck DONALD&quot;</p><p><b>per
former</b>: <a href="organization-example-lab.html">Organization/1832473e-2
fe0-452d-abe9-3cdb9879522f: Acme Laboratory, Inc</a> &quot;Clinical Lab&quo
t;</p><p><b>value</b>: 7.4 mmol/L</p><h3>ReferenceRanges</h3><table class
="grid"><tr><td>-</td><td><b>Low</b></td><td><b>High</b></td></tr><tr><td>*
</td><td>3.0 mmol/L</td><td>7.7 mmol/L</td></tr></table></div>'},
'status': 'final',
'code': {'coding': [{'system': 'http://loinc.org',
'code': '14749-6',
'display': 'Glucose [Moles/volume] in Serum or Plasma'}],
'system': 'http://snomed.info/sct', 'code': '36048009'}}]},
'subject': {'reference': 'Patient/pat2'},
'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
'display': 'Acme Laboratory, Inc'}],
'valueQuantity': {'value': 7.4, 'unit': 'mmol/L'},
'referenceRange': [{'low': {'value': 3, 'unit': 'mmol/L'},
'high': {'value': 7.7, 'unit': 'mmol/L'}}]},
'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801539',
'resource': {'resourceType': 'Observation',
'id': '44801539',
'meta': {'versionId': '1',
'lastUpdated': '2024-07-11T15:44:08.135+00:00',
'source': '#tYPNjWj0MUUsuCAKv'},
'text': {'status': 'generated',
'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narr
ative: Observation</b><a name="o7"> </a></p><div style="display: inline-blo
ck; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px solid
#8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0px
">Resource Observation &quot;o7&quot; </p></div><p><b>status</b>: final</p>
```

```
<p><b>code</b>: Urea nitrogen [Moles/volume] in Serum or Plasma <span style
="background: LightGoldenRodYellow; margin: 4px; border: 1px solid khaki">
(<a href="https://loinc.org/">LOINC</a>#14937-7; <a href="https://browser.i
htsdotools.org/">SNOMED CT</a>#273967009)</span></p><p><b>subject</b>: <a h
ref="patient-example-b.html">Patient/pat2</a> &quot;Duck DONALD&quot;</p><p>
<b>performer</b>: <a href="organization-example-lab.html">Organization/183
2473e-2fe0-452d-abe9-3cdb9879522f: Acme Laboratory, Inc</a> &quot;Clinical
Lab&quot;</p><p><b>value</b>: 4.7 mmol/L</p><h3>ReferenceRanges</h3><table
class="grid"><tr><td>-</td><td><b>Low</b></td><td><b>High</b></td></tr><tr>
<td>*</td><td>2.0 mmol/L</td><td>7.0 mmol/L</td></tr></table></div>}',
  'status': 'final',
  'code': {'coding': [{'system': 'http://loinc.org',
    'code': '14937-7',
    'display': 'Urea nitrogen [Moles/volume] in Serum or Plasma'},
    {'system': 'http://snomed.info/sct', 'code': '273967009'}]},
  'subject': {'reference': 'Patient/pat2'},
  'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
    'display': 'Acme Laboratory, Inc'}],
  'valueQuantity': {'value': 4.7, 'unit': 'mmol/L'},
  'referenceRange': [{'low': {'value': 2, 'unit': 'mmol/L'},
    'high': {'value': 7, 'unit': 'mmol/L'}}],
  'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801540',
  'resource': {'resourceType': 'Observation',
    'id': '44801540',
    'meta': {'versionId': '1',
      'lastUpdated': '2024-07-11T15:44:08.135+00:00',
      'source': '#tYPNjWj0MUuCAKv'},
    'text': {'status': 'generated',
      'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narr
ative: Observation</b><a name="o8"> </a></p><div style="display: inline-blo
ck; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px solid
#8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0px
">Resource Observation &quot;o8&quot;</p></div><p><b>status</b>: final</p>
<p><b>code</b>: Creatinine [Moles/volume] in Serum or Plasma <span style="b
ackground: LightGoldenRodYellow; margin: 4px; border: 1px solid khaki"> (<a
href="https://loinc.org/">LOINC</a>#14682-9; <a href="https://browser.ihtsd
otools.org/">SNOMED CT</a>#70901006)</span></p><p><b>subject</b>: <a href
="patient-example-b.html">Patient/pat2</a> &quot;Duck DONALD&quot;</p><p><b>
performer</b>: <a href="organization-example-lab.html">Organization/183247
3e-2fe0-452d-abe9-3cdb9879522f: Acme Laboratory, Inc</a> &quot;Clinical Lab
&quot;</p><p><b>value</b>: 0.09 mmol/L</p><h3>ReferenceRanges</h3><table cl
ass="grid"><tr><td>-</td><td><b>Low</b></td><td><b>High</b></td></tr><tr><t
d>*</td><td>0.04 mmol/L</td><td>0.11 mmol/L</td></tr></table></div>}',
    'status': 'final',
    'code': {'coding': [{'system': 'http://loinc.org',
      'code': '14682-9',
      'display': 'Creatinine [Moles/volume] in Serum or Plasma'},
        {'system': 'http://snomed.info/sct', 'code': '70901006'}]},
    'subject': {'reference': 'Patient/pat2'},
    'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
      'display': 'Acme Laboratory, Inc'}],
    'valueQuantity': {'value': 0.09, 'unit': 'mmol/L'},
    'referenceRange': [{'low': {'value': 0.04, 'unit': 'mmol/L'},
```

```

    'high': {'value': 0.11, 'unit': 'mmol/L'}}}],
    'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801541',
      'resource': {'resourceType': 'Observation',
        'id': '44801541',
        'meta': {'versionId': '1',
          'lastUpdated': '2024-07-11T15:44:08.135+00:00',
          'source': '#tYPNjWj0MUuCAKv'},
        'text': {'status': 'generated',
          'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narr
ative: Observation</b><a name="o17"> </a></p><div style="display: inline-bl
ock; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px soli
d #8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0
px">Resource Observation &quot;o17&quot; </p></div><p><b>status</b>: final
</p><p><b>code</b>: Calcium [Moles/volume] in Serum or Plasma <span style="
background: LightGoldenRodYellow; margin: 4px; border: 1px solid khaki"> (<
a href="https://loinc.org/">LOINC</a>#2000-8; <a href="https://browser.ihts
dotools.org/">SNOMED CT</a>#71878006)</span></p><p><b>subject</b>: <a href
="patient-example-b.html">Patient/pat2</a> &quot;Duck DONALD&quot;</p><p><b>per
former</b>: <a href="organization-example-lab.html">Organization/183247
3e-2fe0-452d-abe9-3cdb9879522f: Acme Laboratory, Inc</a> &quot;Clinical Lab
&quot;</p><p><b>value</b>: 2.38 mmol/L</p><h3>ReferenceRanges</h3><table cl
ass="grid"><tr><td>-</td><td><b>Low</b></td><td><b>High</b></td></tr><tr><t
d>*</td><td>2.25 mmol/L</td><td>2.65 mmol/L</td></tr></table></div>'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '2000-8',
          'display': 'Calcium [Moles/volume] in Serum or Plasma'},
            {'system': 'http://snomed.info/sct', 'code': '71878006'}]},
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}],
        'valueQuantity': {'value': 2.38, 'unit': 'mmol/L'},
        'referenceRange': [{'low': {'value': 2.25, 'unit': 'mmol/L'},
          'high': {'value': 2.65, 'unit': 'mmol/L'}}]},
        'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801542',
      'resource': {'resourceType': 'Observation',
        'id': '44801542',
        'meta': {'versionId': '1',
          'lastUpdated': '2024-07-11T15:44:08.135+00:00',
          'source': '#tYPNjWj0MUuCAKv'},
        'text': {'status': 'generated',
          'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narr
ative: Observation</b><a name="o1"> </a></p><div style="display: inline-blo
ck; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px solid
#8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0px
">Resource Observation &quot;o1&quot; </p></div><p><b>status</b>: final</p>
<p><b>code</b>: Sodium [Moles/volume] in Serum or Plasma <span style="backg
round: LightGoldenRodYellow; margin: 4px; border: 1px solid khaki"> (<a hre
f="https://loinc.org/">LOINC</a>#2951-2; <a href="https://browser.ihtsdotoo
ls.org/">SNOMED CT</a>#104934005)</span></p><p><b>subject</b>: <a href="pat
ient-example-b.html">Patient/pat2</a> &quot;Duck DONALD&quot;</p><p><b>perf
ormer</b>: <a href="organization-example-lab.html">Organization/1832473e-2f
e0-452d-abe9-3cdb9879522f: Acme Laboratory, Inc</a> &quot;Clinical Lab&quo

```

```
t;</p><p><b>value</b>: 140 mmol/L</p><h3>ReferenceRanges</h3><table class
="grid"><tr><td></td><td><b>Low</b></td><td><b>High</b></td></tr><tr><td>*
</td><td>137 mmol/L</td><td>147 mmol/L</td></tr></table></div>'},
  'status': 'final',
  'code': {'coding': [{'system': 'http://loinc.org',
    'code': '2951-2',
    'display': 'Sodium [Moles/volume] in Serum or Plasma'},
    {'system': 'http://snomed.info/sct', 'code': '104934005'}]},
  'subject': {'reference': 'Patient/pat2'},
  'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
    'display': 'Acme Laboratory, Inc'}]},
  'valueQuantity': {'value': 140, 'unit': 'mmol/L'},
  'referenceRange': [{'low': {'value': 137, 'unit': 'mmol/L'},
    'high': {'value': 147, 'unit': 'mmol/L'}}]},
  'search': {'mode': 'match'}},
  {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801543',
  'resource': {'resourceType': 'Observation',
    'id': '44801543',
    'meta': {'versionId': '1',
      'lastUpdated': '2024-07-11T15:44:08.135+00:00',
      'source': '#tYPNjWj0MUuCAKv'},
    'text': {'status': 'generated',
      'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narr
ative: Observation</b><a name="o2"> </a></p><div style="display: inline-blo
ck; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px solid
#8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0px
">Resource Observation &quot;o2&quot; </p></div><p><b>status</b>: final</p>
<p><b>code</b>: Potassium [Moles/volume] in Serum or Plasma <span style="ba
ckground: LightGoldenRodYellow; margin: 4px; border: 1px solid khaki"> (<a
href="https://loinc.org/">LOINC</a>#2823-3; <a href="https://browser.ihtsdo
tools.org/">SNOMED CT</a>#59573005)</span></p><p><b>subject</b>: <a href="p
atient-example-b.html">Patient/pat2</a> &quot;Duck DONALD&quot;</p><p><b>pe
rformer</b>: <a href="organization-example-lab.html">Organization/1832473
e-2fe0-452d-abe9-3cdb9879522f: Acme Laboratory, Inc</a> &quot;Clinical Lab&
quot;</p><p><b>value</b>: 4.2 mmol/L</p><h3>ReferenceRanges</h3><table clas
s="grid"><tr><td></td><td><b>Low</b></td><td><b>High</b></td></tr><tr><td>
*</td><td>3.5 mmol/L</td><td>5.0 mmol/L</td></tr></table></div>'},
    'status': 'final',
    'code': {'coding': [{'system': 'http://loinc.org',
      'code': '2823-3',
      'display': 'Potassium [Moles/volume] in Serum or Plasma'},
      {'system': 'http://snomed.info/sct', 'code': '59573005'}]},
    'subject': {'reference': 'Patient/pat2'},
    'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
      'display': 'Acme Laboratory, Inc'}]},
    'valueQuantity': {'value': 4.2, 'unit': 'mmol/L'},
    'referenceRange': [{'low': {'value': 3.5, 'unit': 'mmol/L'},
      'high': {'value': 5, 'unit': 'mmol/L'}}]},
    'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801544',
    'resource': {'resourceType': 'Observation',
      'id': '44801544',
      'meta': {'versionId': '1',
        'lastUpdated': '2024-07-11T15:44:08.135+00:00',
```



```
'source': '#tYPNjWj0MUsuCAKv'},
'text': {'status': 'generated',
'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narrative: Observation</b><a name="o3"> </a></p><div style="display: inline-block; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px solid #8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0px">Resource Observation &quot;o3&quot; </p></div><p><b>status</b>: final</p><p><b>code</b>: Chloride [Moles/volume] in Serum or Plasma <span style="background: LightGoldenRodYellow; margin: 4px; border: 1px solid khaki"> (<a href="https://loinc.org/">LOINC</a>#2075-0; <a href="https://browser.ihtsdotools.org/">SNOMED CT</a>#46511006)</span></p><p><b>subject</b>: <a href="patient-example-b.html">Patient/pat2</a> &quot;Duck DONALD&quot; </p><p><b>performer</b>: <a href="organization-example-lab.html">Organization/1832473e-2fe0-452d-abe9-3cdb9879522f: Acme Laboratory, Inc</a> &quot;Clinical Lab&quot; </p><p><b>value</b>: 105 mmol/L</p><h3>ReferenceRanges</h3><table class="grid"><tr><td>-</td><td><b>Low</b></td><td><b>High</b></td></tr><tr><td>*</td><td>96 mmol/L</td><td>109 mmol/L</td></tr></table></div>',
'status': 'final',
'code': {'coding': [{'system': 'http://loinc.org',
'code': '2075-0',
'display': 'Chloride [Moles/volume] in Serum or Plasma'},
{'system': 'http://snomed.info/sct', 'code': '46511006'}]},
'subject': {'reference': 'Patient/pat2'},
'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9879522f',
'display': 'Acme Laboratory, Inc'}],
'valueQuantity': {'value': 105, 'unit': 'mmol/L'},
'referenceRange': [{'low': {'value': 96, 'unit': 'mmol/L'},
'high': {'value': 109, 'unit': 'mmol/L'}}],
'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801545',
'resource': {'resourceType': 'Observation',
'id': '44801545',
'meta': {'versionId': '1',
'lastUpdated': '2024-07-11T15:44:08.135+00:00',
'source': '#tYPNjWj0MUsuCAKv'},
'text': {'status': 'generated',
'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narrative: Observation</b><a name="o4"> </a></p><div style="display: inline-block; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px solid #8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0px">Resource Observation &quot;o4&quot; </p></div><p><b>status</b>: final</p><p><b>code</b>: Bicarbonate [Moles/volume] in Serum <span style="background: LightGoldenRodYellow; margin: 4px; border: 1px solid khaki"> (<a href="https://loinc.org/">LOINC</a>#1963-8; <a href="https://browser.ihtsdotools.org/">SNOMED CT</a>#88645003)</span></p><p><b>subject</b>: <a href="patient-example-b.html">Patient/pat2</a> &quot;Duck DONALD&quot; </p><p><b>performer</b>: <a href="organization-example-lab.html">Organization/1832473e-2fe0-452d-abe9-3cdb9879522f: Acme Laboratory, Inc</a> &quot;Clinical Lab&quot; </p><p><b>value</b>: 26 mmol/L</p><h3>ReferenceRanges</h3><table class="grid"><tr><td>-</td><td><b>Low</b></td><td><b>High</b></td></tr><tr><td>*</td><td>25 mmol/L</td><td>33 mmol/L</td></tr></table></div>',
'status': 'final',
'code': {'coding': [{'system': 'http://loinc.org',
'code': '1963-8',
'display': 'Bicarbonate [Moles/volume] in Serum'},
```

```

      {'system': 'http://snomed.info/sct', 'code': '88645003'}}]],
      'subject': {'reference': 'Patient/pat2'},
      'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
        'display': 'Acme Laboratory, Inc'}],
      'valueQuantity': {'value': 26, 'unit': 'mmol/L'},
      'referenceRange': [{'low': {'value': 25, 'unit': 'mmol/L'},
        'high': {'value': 33, 'unit': 'mmol/L'}}]],
      'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801546',
      'resource': {'resourceType': 'Observation',
        'id': '44801546',
        'meta': {'versionId': '1',
          'lastUpdated': '2024-07-11T15:51:37.783+00:00',
          'source': '#kZ3KAReIzyBIPCpk'},
        'text': {'status': 'generated',
          'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narr
ative: Observation</b><a name="o14"> </a></p><div style="display: inline-bl
ock; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px soli
d #8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0
px">Resource Observation &quot;o14&quot; </p></div><p><b>status</b>: final
</p><p><b>code</b>: Alanine aminotransferase [Enzymatic activity/volume] in
Serum or Plasma <span style="background: LightGoldenRodYellow; margin: 4px;
border: 1px solid khaki"> (<a href="https://loinc.org/">LOINC</a>#1742-6; <
a href="https://browser.ihtsdotools.org/">SNOMED CT</a>#34608000)</span></p>
<p><b>subject</b>: <a href="patient-example-b.html">Patient/pat2</a> &quot;
;Duck DONALD&quot;</p><p><b>performer</b>: <a href="organization-example-la
b.html">Organization/1832473e-2fe0-452d-abe9-3cdb9879522f: Acme Laboratory,
Inc</a> &quot;Clinical Lab&quot;</p><p><b>value</b>: 54 U/L</p><p><b>interp
retation</b>: High <span style="background: LightGoldenRodYellow; margin: 4
px; border: 1px solid khaki"> (<a href="http://terminology.hl7.org/5.1.0/Co
deSystem-v3-ObservationInterpretation.html">ObservationInterpretation</a>#
H)</span></p><h3>ReferenceRanges</h3><table class="grid"><tr><td>-</td><td>
<b>Low</b></td><td><b>High</b></td></tr><tr><td>*</td><td>0 U/L</td><td>45
U/L</td></tr></table></div>'},
        'status': 'final',
        'code': {'coding': [{'system': 'http://loinc.org',
          'code': '1742-6',
          'display': 'Alanine aminotransferase [Enzymatic activity/volume] in
Serum or Plasma'}],
          {'system': 'http://snomed.info/sct', 'code': '34608000'}}]],
        'subject': {'reference': 'Patient/pat2'},
        'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
          'display': 'Acme Laboratory, Inc'}],
        'valueQuantity': {'value': 54, 'unit': 'U/L'},
        'interpretation': [{'coding': [{'system': 'http://terminology.hl7.org/C
odeSystem/v3-ObservationInterpretation',
          'code': 'H'}}]],
        'referenceRange': [{'low': {'value': 0, 'unit': 'U/L'},
          'high': {'value': 45, 'unit': 'U/L'}}]],
        'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801547',
      'resource': {'resourceType': 'Observation',
        'id': '44801547',
        'meta': {'versionId': '1',

```

```

    'lastUpdated': '2024-07-11T15:51:37.783+00:00',
    'source': '#kZ3KAReIzyBIPCpk'},
    'text': {'status': 'generated',
      'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narr
ative: Observation</b><a name="o15"> </a></p><div style="display: inline-bl
ock; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px soli
d #8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0
px">Resource Observation &quot;o15&quot; </p></div><p><b>status</b>: final
</p><p><b>code</b>: Aspartate aminotransferase [Enzymatic activity/volume]
in Serum or Plasma <span style="background: LightGoldenRodYellow; margin: 4
px; border: 1px solid khaki"> (<a href="https://loinc.org/">LOINC</a>#192
0-8; <a href="https://browser.ihtsdotools.org/">SNOMED CT</a>#45896001)</sp
an></p><p><b>subject</b>: <a href="patient-example-b.html">Patient/pat2</a>
&quot;Duck DONALD&quot;</p><p><b>performer</b>: <a href="organization-examp
le-lab.html">Organization/1832473e-2fe0-452d-abe9-3cdb9879522f: Acme Labora
tory, Inc</a> &quot;Clinical Lab&quot;</p><p><b>value</b>: 30 U/L</p><h3>Re
ferenceRanges</h3><table class="grid"><tr><td>-</td><td><b>Low</b></td><td>
<b>High</b></td></tr><tr><td>*</td><td>0 U/L</td><td>41 U/L</td></tr></tabl
e></div>'},
    'status': 'final',
    'code': {'coding': [{'system': 'http://loinc.org',
      'code': '1920-8',
      'display': 'Aspartate aminotransferase [Enzymatic activity/volume] i
n Serum or Plasma'}],
      'system': 'http://snomed.info/sct', 'code': '45896001'}]},
    'subject': {'reference': 'Patient/pat2'},
    'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
      'display': 'Acme Laboratory, Inc'}],
    'valueQuantity': {'value': 30, 'unit': 'U/L'},
    'referenceRange': [{'low': {'value': 0, 'unit': 'U/L'},
      'high': {'value': 41, 'unit': 'U/L'}}],
    'search': {'mode': 'match'}},
    {'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801548',
      'resource': {'resourceType': 'Observation',
        'id': '44801548',
        'meta': {'versionId': '1',
          'lastUpdated': '2024-07-11T15:51:37.783+00:00',
          'source': '#kZ3KAReIzyBIPCpk'},
          'text': {'status': 'generated',
            'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narr
ative: Observation</b><a name="o6"> </a></p><div style="display: inline-blo
ck; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px solid
#8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0px
">Resource Observation &quot;o6&quot; </p></div><p><b>status</b>: final</p>
<p><b>code</b>: Glucose [Moles/volume] in Serum or Plasma <span style="back
ground: LightGoldenRodYellow; margin: 4px; border: 1px solid khaki"> (<a hr
ef="https://loinc.org/">LOINC</a>#14749-6; <a href="https://browser.ihtsdot
ools.org/">SNOMED CT</a>#36048009)</span></p><p><b>subject</b>: <a href="pa
tient-example-b.html">Patient/pat2</a> &quot;Duck DONALD&quot;</p><p><b>per
former</b>: <a href="organization-example-lab.html">Organization/1832473e-2
fe0-452d-abe9-3cdb9879522f: Acme Laboratory, Inc</a> &quot;Clinical Lab&quo
t;</p><p><b>value</b>: 7.4 mmol/L</p><h3>ReferenceRanges</h3><table class
="grid"><tr><td>-</td><td><b>Low</b></td><td><b>High</b></td></tr><tr><td>*</td>
<td>3.0 mmol/L</td><td>7.7 mmol/L</td></tr></table></div>'},
          'status': 'final',

```

```

'code': {'coding': [{'system': 'http://loinc.org',
  'code': '14749-6',
  'display': 'Glucose [Moles/volume] in Serum or Plasma'}],
  {'system': 'http://snomed.info/sct', 'code': '36048009'}}],
'subject': {'reference': 'Patient/pat2'},
'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
  'display': 'Acme Laboratory, Inc'}],
'valueQuantity': {'value': 7.4, 'unit': 'mmol/L'},
'referenceRange': [{'low': {'value': 3, 'unit': 'mmol/L'},
  'high': {'value': 7.7, 'unit': 'mmol/L'}}],
'search': {'mode': 'match'}},
{'fullUrl': 'https://hapi.fhir.org/baseR4/Observation/44801549',
'resource': {'resourceType': 'Observation',
'id': '44801549',
'meta': {'versionId': '1',
'lastUpdated': '2024-07-11T15:51:37.783+00:00',
'source': '#kZ3KAReIzyBIPCpk'},
'text': {'status': 'generated',
'div': '<div xmlns="http://www.w3.org/1999/xhtml"><p><b>Generated Narr
ative: Observation</b><a name="o7"> </a></p><div style="display: inline-blo
ck; background-color: #d9e0e7; padding: 6px; margin: 4px; border: 1px solid
#8da1b4; border-radius: 5px; line-height: 60%"><p style="margin-bottom: 0px
">Resource Observation &quot;o7&quot; </p></div><p><b>status</b>: final</p>
<p><b>code</b>: Urea nitrogen [Moles/volume] in Serum or Plasma <span style
="background: LightGoldenRodYellow; margin: 4px; border: 1px solid khaki">
(<a href="https://loinc.org/">LOINC</a>#14937-7; <a href="https://browser.i
htsdotools.org/">SNOMED CT</a>#273967009)</span></p><p><b>subject</b>: <a h
ref="patient-example-b.html">Patient/pat2</a> &quot;Duck DONALD&quot;</p><p>
<b>performer</b>: <a href="organization-example-lab.html">Organization/183
2473e-2fe0-452d-abe9-3cdb9879522f: Acme Laboratory, Inc</a> &quot;Clinical
Lab&quot;</p><p><b>value</b>: 4.7 mmol/L</p><h3>ReferenceRanges</h3><table
class="grid"><tr><td>-</td><td><b>Low</b></td><td><b>High</b></td></tr><tr>
<td>*</td><td>2.0 mmol/L</td><td>7.0 mmol/L</td></tr></table></div>'},
'status': 'final',
'code': {'coding': [{'system': 'http://loinc.org',
  'code': '14937-7',
  'display': 'Urea nitrogen [Moles/volume] in Serum or Plasma'}],
  {'system': 'http://snomed.info/sct', 'code': '273967009'}}],
'subject': {'reference': 'Patient/pat2'},
'performer': [{'reference': 'Organization/1832473e-2fe0-452d-abe9-3cdb9
879522f',
  'display': 'Acme Laboratory, Inc'}],
'valueQuantity': {'value': 4.7, 'unit': 'mmol/L'},
'referenceRange': [{'low': {'value': 2, 'unit': 'mmol/L'},
  'high': {'value': 7, 'unit': 'mmol/L'}}],
'search': {'mode': 'match'}}}]

```

In [5]: *# let's write some code to unpack this. the results come back as a bundle and
can use our json viewer to see that part of this is a list of resources
in each resource we can find the lab names and values!*

```
observations = []
```

```
for entry in obs.get("entry", []):
    resource = entry["resource"]
```

```

if resource.get("code", {}).get("coding"):
    code = resource.get("code", {}).get("coding")[0].get("display")
    value = resource.get("valueQuantity", {}).get("value")
    unit = resource.get("valueQuantity", {}).get("unit")
    time = resource.get("effectiveDateTime")
else:
    print("Skipping entry due to missing 'coding' information.")

if value:
    observations.append({"code": code, "value": value, "unit": unit, "ti

df = pd.DataFrame(observations)
print(df.head())

```

Skipping entry due to missing 'coding' information.

		code	value	unit	time
0	Cholesterol [Moles/volume] in Serum or Plasma		6.3	mmol/L	None
1	Hemoglobin [Mass/volume] in Blood		176.0	g/L	None
2	Erythrocytes [# /volume] in Blood by Automated ...		5.9	x10*12/L	None
3	Hematocrit [Volume Fraction] of Blood by Autom...		55.0	%	None
4	Erythrocyte mean corpuscular volume [Entitic v...		99.0	fL	None

In [6]: *# as you can see this is a random collection of labs... for our app we need
detail.*

```

# let's create our own patient now so we have control over the data
# we are allowed to do this because it is a public endpoint
# this patient will be purged at some interval

```

In [7]: *# this makes a fake new patient*
if everything goes well, we'll get a 201 response and an ID

```

new_patients = [
    {
        "resourceType": "Patient",
        "name": [{"use": "official", "family": "Doe", "given": ["Jane"]}],
        "gender": "female",
        "birthDate": "1980-05-01"
    },
    {
        'resourceType': 'Patient',
        'name': [{'use': 'official', 'family': 'Smith', 'given': ['Martina']
        'gender': 'female',
        'birthDate': '2003-12-05'
    }
]

for patient in new_patients:
    resp = requests.post(f"{FHIR_BASE}/Patient", json=patient)
    print(resp.status_code, resp.json())
    patient_id = resp.json()["id"]
    print("New Patient ID:", patient_id)

```

```
201 {'resourceType': 'Patient', 'id': '50166462', 'meta': {'versionId': '1',
'lastUpdated': '2025-10-14T19:32:26.538+00:00', 'source': '#SLK7GQSXIsyva99
D'}, 'text': {'status': 'generated', 'div': '<div xmlns="http://www.w3.org/1
999/xhtml"><div class="hapiHeaderText">Jane <b>DOE </b></div><table class="h
apiPropertyTable"><tbody><tr><td>Date of birth</td><td><span>01 May 1980</sp
an></td></tr></tbody></table></div>'}, 'name': [{'use': 'official', 'family'
: 'Doe', 'given': ['Jane']}], 'gender': 'female', 'birthDate': '1980-05-01'}
New Patient ID: 50166462
```

```
201 {'resourceType': 'Patient', 'id': '50166463', 'meta': {'versionId': '1',
'lastUpdated': '2025-10-14T19:32:27.539+00:00', 'source': '#muB8mQMNTyb8xEU
K'}, 'text': {'status': 'generated', 'div': '<div xmlns="http://www.w3.org/1
999/xhtml"><div class="hapiHeaderText">Martina <b>SMITH </b></div><table cla
ss="hapiPropertyTable"><tbody><tr><td>Date of birth</td><td><span>05 Decembe
r 2003</span></td></tr></tbody></table></div>'}, 'name': [{'use': 'official'
, 'family': 'Smith', 'given': ['Martina']}], 'gender': 'female', 'birthDate'
: '2003-12-05'}
New Patient ID: 50166463
```

```
In [8]: patient_id = resp.json()["id"]
print("New Patient ID:", patient_id)
```

New Patient ID: 50166463

```
In [9]: # let's add a bunch of random data, one set of observations per day

from random import randint, gauss
from datetime import datetime, timedelta

base_date = datetime(2025, 8, 1)

observations = []
for i in range(5): # 5 days
    dt = (base_date + timedelta(days=i)).isoformat() + "Z"

    observations += [
        {
            "resourceType": "Observation",
            "status": "final",
            "category": [{"coding": [{"system": "http://terminology.hl7.org/
                                "code": "vital-signs"}]}],
            "code": {"text": "Systolic Blood Pressure"},
            "subject": {"reference": f"Patient/{patient_id}"},
            "effectiveDateTime": dt,
            "valueQuantity": {"value": randint(110, 135), "unit": "mmHg"}
        },
        {
            "resourceType": "Observation",
            "status": "final",
            "category": [{"coding": [{"system": "http://terminology.hl7.org/
                                "code": "vital-signs"}]}],
            "code": {"text": "Diastolic Blood Pressure"},
            "subject": {"reference": f"Patient/{patient_id}"},
            "effectiveDateTime": dt,
            "valueQuantity": {"value": randint(70, 90), "unit": "mmHg"}
        },
        {
            "resourceType": "Observation",
```

```

        "status": "final",
        "category": [{"coding": [{"system": "http://terminology.hl7.org/
                                "code": "vital-signs"}]}],
        "code": {"text": "Heart rate"},
        "subject": {"reference": f"Patient/{patient_id}"},
        "effectiveDateTime": dt,
        "valueQuantity": {"value": randint(65, 90), "unit": "beats/min"}
    },
    {
        "resourceType": "Observation",
        "status": "final",
        "category": [{"coding": [{"system": "http://terminology.hl7.org/
                                "code": "laboratory"}]}],
        "code": {"text": "Glucose"},
        "subject": {"reference": f"Patient/{patient_id}"},
        "effectiveDateTime": dt,
        "valueQuantity": {"value": round(gauss(95, 10), 1), "unit": "mg/
    }
]

# POST them
for obs in observations:
    r = requests.post(f"{FHIR_BASE}/Observation", json=obs, params={"patient": patient_id})
    print("Obs:", r.status_code, r.json().get("id"))

```

```

Obs: 201 50166465
Obs: 201 50166466
Obs: 201 50166467
Obs: 201 50166468
Obs: 201 50166469
Obs: 201 50166471
Obs: 201 50166472
Obs: 201 50166473
Obs: 201 50166474
Obs: 201 50166475
Obs: 201 50166477
Obs: 201 50166478
Obs: 201 50166479
Obs: 201 50166480
Obs: 201 50166481
Obs: 201 50166483
Obs: 201 50166484
Obs: 201 50166485
Obs: 201 50166486
Obs: 201 50166487

```

```

In [10]: # Query observations
obs_resp = requests.get(f"{FHIR_BASE}/Observation", params={"patient": patient_id})
print(len(obs_resp.json().get("entry", [])), "observations found")
print(patient_id)

```

```

20 observations found
50166463

```

```

In [11]: observations = []
for entry in obs_resp.json().get("entry", []):
    resource = entry["resource"]

```

```

code = resource.get("code", {}).get("text")
value = resource.get("valueQuantity", {}).get("value")
unit = resource.get("valueQuantity", {}).get("unit")
time = resource.get("effectiveDateTime")
if value:
    observations.append({"code": code, "value": value, "unit": unit, "ti

df = pd.DataFrame(observations)
print(df)

```

		code	value	unit	time
0	Systolic Blood Pressure	131.0	mmHg	2025-08-01T00:00:00Z	
1	Diastolic Blood Pressure	72.0	mmHg	2025-08-01T00:00:00Z	
2	Heart rate	75.0	beats/min	2025-08-01T00:00:00Z	
3	Glucose	100.6	mg/dL	2025-08-01T00:00:00Z	
4	Systolic Blood Pressure	110.0	mmHg	2025-08-02T00:00:00Z	
5	Diastolic Blood Pressure	79.0	mmHg	2025-08-02T00:00:00Z	
6	Heart rate	78.0	beats/min	2025-08-02T00:00:00Z	
7	Glucose	99.5	mg/dL	2025-08-02T00:00:00Z	
8	Systolic Blood Pressure	122.0	mmHg	2025-08-03T00:00:00Z	
9	Diastolic Blood Pressure	84.0	mmHg	2025-08-03T00:00:00Z	
10	Heart rate	71.0	beats/min	2025-08-03T00:00:00Z	
11	Glucose	84.3	mg/dL	2025-08-03T00:00:00Z	
12	Systolic Blood Pressure	111.0	mmHg	2025-08-04T00:00:00Z	
13	Diastolic Blood Pressure	88.0	mmHg	2025-08-04T00:00:00Z	
14	Heart rate	89.0	beats/min	2025-08-04T00:00:00Z	
15	Glucose	110.0	mg/dL	2025-08-04T00:00:00Z	
16	Systolic Blood Pressure	116.0	mmHg	2025-08-05T00:00:00Z	
17	Diastolic Blood Pressure	90.0	mmHg	2025-08-05T00:00:00Z	
18	Heart rate	83.0	beats/min	2025-08-05T00:00:00Z	
19	Glucose	89.0	mg/dL	2025-08-05T00:00:00Z	

```

In [12]: medications = [
    {
        "resourceType": "MedicationStatement",
        "status": "active",
        "medicationCodeableConcept": {"text": "Lisinopril 10mg tablet"},
        "subject": {"reference": f"Patient/{patient_id}"},
        "effectiveDateTime": "2025-09-10T11:00:00Z"
    },
    {
        "resourceType": "MedicationStatement",
        "status": "active",
        "medicationCodeableConcept": {"text": "Lisinopril 10mg tablet"},
        "subject": {"reference": f"Patient/{patient_id}"},
        "effectiveDateTime": "2025-09-11T11:30:00Z"
    },
    {
        "resourceType": "MedicationStatement",
        "status": "active",
        "medicationCodeableConcept": {"text": "Lisinopril 20mg tablet"},
        "subject": {"reference": f"Patient/{patient_id}"},
        "effectiveDateTime": "2025-09-12T12:00:00Z"
    }
]

```



```

for med in medications:
    r = requests.post(f"{FHIR_BASE}/MedicationStatement", json=med)
    print("Med:", r.status_code, r.json().get("id"))

```

Med: 201 50166489

Med: 201 50166490

Med: 201 50166491

In [13]: patient_id

Out[13]: '50166463'

In [14]: *# Now that we have our fake patient set up, let's go back to the goal and put their data*

Observations

```
obs_resp = requests.get(f"{FHIR_BASE}/Observation", params={"patient": patient_id})
```

```
print("----- OBSERVATIONS -----")
```

```
print(len(obs_resp.json().get("entry", [])), "observations")
```

```
observations = []
```

```
for entry in obs_resp.json().get("entry", []):
```

```
    resource = entry["resource"]
```

```
    code = resource.get("code", {}).get("text")
```

```
    value = resource.get("valueQuantity", {}).get("value")
```

```
    unit = resource.get("valueQuantity", {}).get("unit")
```

```
    time = resource.get("effectiveDateTime")
```

```
    if value:
```

```
        observations.append({"code": code, "value": value, "unit": unit, "time": time})
```

```
df = pd.DataFrame(observations)
```

```
print(df)
```

Medications

```
print("\n----- MEDICATIONS -----")
```

```
med_resp = requests.get(f"{FHIR_BASE}/MedicationStatement", params={"patient": patient_id})
```

```
for e in med_resp.json().get("entry", []):
```

```
    print("Medication:", e["resource"].get("medicationCodeableConcept", {}).
```

----- OBSERVATIONS -----

20 observations

		code	value	unit	time
0	Systolic Blood Pressure		131.0	mmHg	2025-08-01T00:00:00Z
1	Diastolic Blood Pressure		72.0	mmHg	2025-08-01T00:00:00Z
2	Heart rate		75.0	beats/min	2025-08-01T00:00:00Z
3	Glucose		100.6	mg/dL	2025-08-01T00:00:00Z
4	Systolic Blood Pressure		110.0	mmHg	2025-08-02T00:00:00Z
5	Diastolic Blood Pressure		79.0	mmHg	2025-08-02T00:00:00Z
6	Heart rate		78.0	beats/min	2025-08-02T00:00:00Z
7	Glucose		99.5	mg/dL	2025-08-02T00:00:00Z
8	Systolic Blood Pressure		122.0	mmHg	2025-08-03T00:00:00Z
9	Diastolic Blood Pressure		84.0	mmHg	2025-08-03T00:00:00Z
10	Heart rate		71.0	beats/min	2025-08-03T00:00:00Z
11	Glucose		84.3	mg/dL	2025-08-03T00:00:00Z
12	Systolic Blood Pressure		111.0	mmHg	2025-08-04T00:00:00Z
13	Diastolic Blood Pressure		88.0	mmHg	2025-08-04T00:00:00Z
14	Heart rate		89.0	beats/min	2025-08-04T00:00:00Z
15	Glucose		110.0	mg/dL	2025-08-04T00:00:00Z
16	Systolic Blood Pressure		116.0	mmHg	2025-08-05T00:00:00Z
17	Diastolic Blood Pressure		90.0	mmHg	2025-08-05T00:00:00Z
18	Heart rate		83.0	beats/min	2025-08-05T00:00:00Z
19	Glucose		89.0	mg/dL	2025-08-05T00:00:00Z

----- MEDICATIONS -----

Medication: Lisinopril 10mg tablet

Medication: Lisinopril 10mg tablet

Medication: Lisinopril 20mg tablet

In [15]: # ok now let's plot!

```

fig, axes = plt.subplots(3, 1, figsize=(10, 8), gridspec_kw={"height_ratios"


# show bp on one plot
bp = df[df["code"].str.contains("Blood Pressure")]
for code in ["Systolic Blood Pressure", "Diastolic Blood Pressure"]:
    subset = bp[bp["code"] == code]
    axes[0].plot(subset["time"], subset["value"], marker="o", label=code)
axes[0].set_ylabel("mmHg")
axes[0].set_title("Blood Pressure Over Time")
axes[0].legend()
axes[0].grid(True)

hr = df[df["code"] == "Heart rate"]
axes[1].plot(hr["time"], hr["value"], marker="o", color="tab:orange", label="")
axes[1].set_ylabel("beats/min")
axes[1].legend()
axes[1].grid(True)

glucose = df[df["code"] == "Glucose"]
axes[2].plot(glucose["time"], glucose["value"], marker="o", color="tab:green")
axes[2].set_ylabel("mg/dL")
axes[2].legend()
axes[2].grid(True)

```

```
plt.show()
```

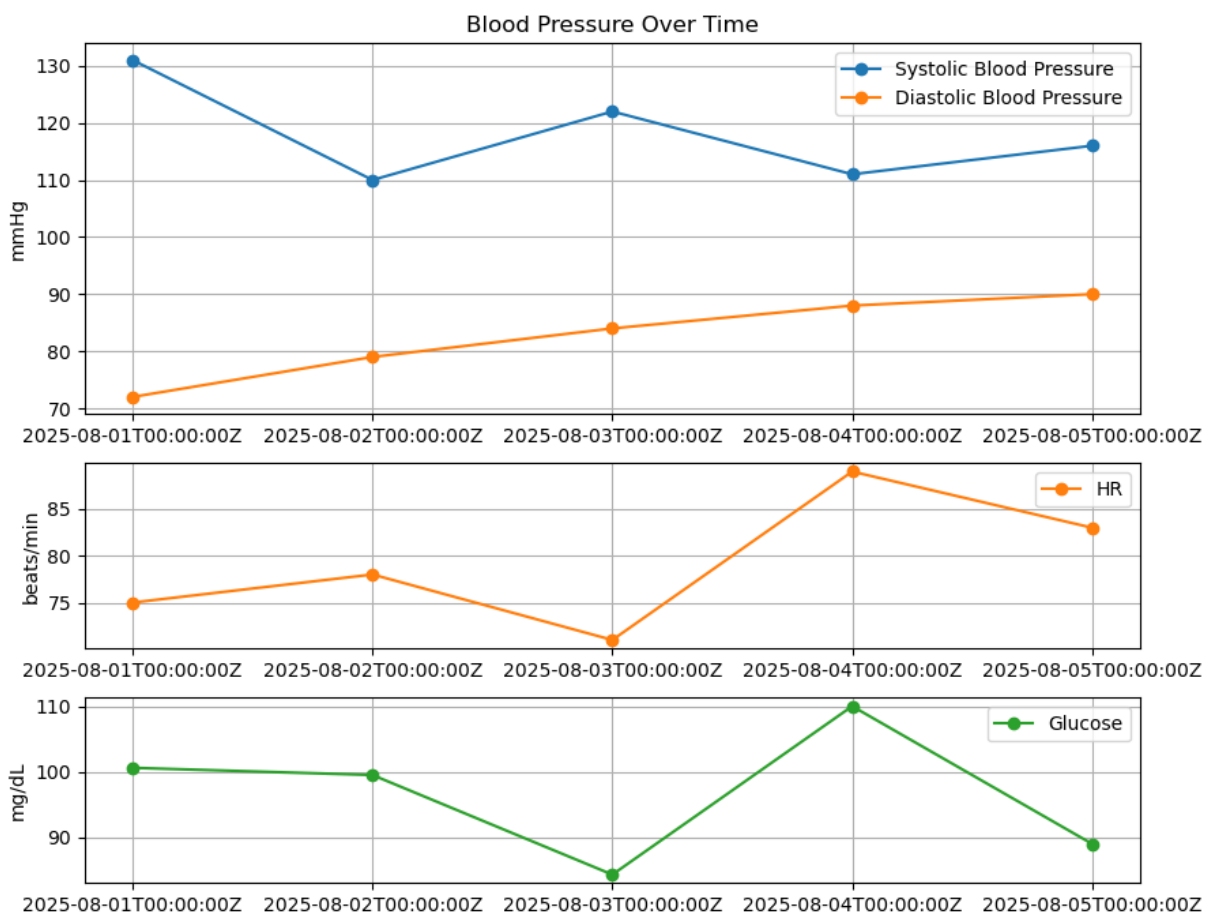

```
In [1]: # if you made it this far congratulations! To extend your knowledge, please
# find a new resource type on the FHIR specification, make some synthetic data
# and retrieve it, showing the data however you like.

# you can browse resources here: https://build.fhir.org/resource-list.html
```

Make synthetic data with immunization

```
In [16]: {
  "coding": [
    {
      "system": "http://hl7.org/fhir/sid/cvx",
      "code": "102",
      "display": "MMR"
    }
  ]
}

Out[16]: {'coding': [{'system': 'http://hl7.org/fhir/sid/cvx',
  'code': '102',
  'display': 'MMR'}]}
```

```
In [17]: {
  "reference": "Patient/49233412"
}
```

```
Out[17]: {'reference': 'Patient/49233412'}
```

```
In [18]: {
  "resourceType": "Immunization",
  "status": "completed",
  "vaccineCode": {
    "coding": [
      {
        "system": "http://hl7.org/fhir/sid/cvx",
        "code": "102",
        "display": "MMR"
      }
    ]
  },
  "patient": {
    "reference": "Patient/49233412"
  },
  "occurrenceDateTime": "2025-08-15T10:00:00Z",
  "performer": [
    {
      "function": {
        "coding": [
          {
            "system": "http://terminology.hl7.org/CodeSystem/v2-0443",
            "code": "AP",
            "display": "Administering Provider"
          }
        ]
      }
    }
  ],
  "actor": {
    "reference": "Practitioner/example-practitioner-id"
  }
},
  "lotNumber": "ABC12345",
  "expirationDate": "2026-08-15"
}
```

```
Out[18]: {'resourceType': 'Immunization',
  'status': 'completed',
  'vaccineCode': {'coding': [{'system': 'http://hl7.org/fhir/sid/cvx',
    'code': '102',
    'display': 'MMR'}]}},
  'patient': {'reference': 'Patient/49233412'},
  'occurrenceDateTime': '2025-08-15T10:00:00Z',
  'performer': [{'function': {'coding': [{'system': 'http://terminology.hl
7.org/CodeSystem/v2-0443',
    'code': 'AP',
    'display': 'Administering Provider'}]}},
  'actor': {'reference': 'Practitioner/example-practitioner-id'}}],
  'lotNumber': 'ABC12345',
  'expirationDate': '2026-08-15'}
```

```
In [19]: import requests
import json
```

```

FHIR_BASE = "http://hapi.fhir.org/baseR4"
patient_id = "49233412"

immunization_data = {
    "resourceType": "Immunization",
    "status": "completed",
    "vaccineCode": {
        "coding": [
            {
                "system": "http://hl7.org/fhir/sid/cvx",
                "code": "102",
                "display": "MMR"
            }
        ],
        "text": "MMR"
    },
    "patient": {
        "reference": f"Patient/{patient_id}"
    },
    "occurrenceDateTime": "2025-08-15T10:00:00Z",
    "lotNumber": "ABC12345"
}

response = requests.post(f"{FHIR_BASE}/Immunization", json=immunization_data)
print(f"Status Code: {response.status_code}")
print(response.json())

```

Status Code: 201

```

{'resourceType': 'Immunization', 'id': '50166507', 'meta': {'versionId': '1',
'lastUpdated': '2025-10-14T19:33:59.477+00:00', 'source': '#ipCD4S7Yi0VF7q0j'},
'status': 'completed', 'vaccineCode': {'coding': [{'system': 'http://hl7.org/fhir/sid/cvx',
'code': '102', 'display': 'MMR'}]}, 'text': 'MMR', 'patient': {'reference': 'Patient/49233412'},
'occurrenceDateTime': '2025-08-15T10:00:00Z', 'lotNumber': 'ABC12345'}

```

In [20]:

```

FHIR_BASE = "http://hapi.fhir.org/baseR4"
immunization_id = "50166507"

response = requests.get(f"{FHIR_BASE}/Immunization/{immunization_id}")
print(response.json())

```

```

{'resourceType': 'Immunization', 'id': '50166507', 'meta': {'versionId': '1',
'lastUpdated': '2025-10-14T19:33:59.477+00:00', 'source': '#ipCD4S7Yi0VF7q0j'},
'status': 'completed', 'vaccineCode': {'coding': [{'system': 'http://hl7.org/fhir/sid/cvx',
'code': '102', 'display': 'MMR'}]}, 'text': 'MMR', 'patient': {'reference': 'Patient/49233412'},
'occurrenceDateTime': '2025-08-15T10:00:00Z', 'lotNumber': 'ABC12345'}

```